


Installation of casters for LiteGait Standard Profile Base

Tools Required:

- 1/2" open ended wrench
- 3/4" open ended wrench
- 5/16" hex key
- Flat head screwdriver

Casters should be installed with two total lock casters on the left leg and two directional lock casters on the right leg.

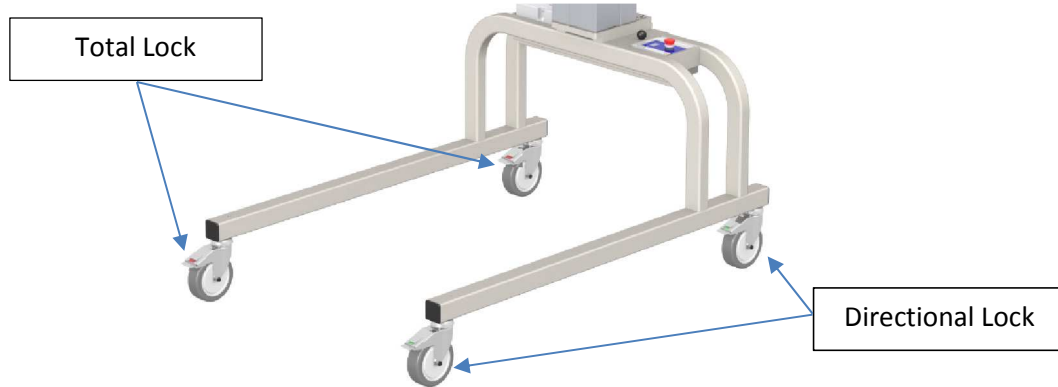

1) Remove all loose accessories, harness, BiSym scale, battery and hand switch before leaning LiteGait to Access casters. Using two people lean LiteGait to access casters.

2) Lock the casters to prevent swiveling of the casters.



3) Use 1/2" Open Ended Wrench to loosen the caster hardware.

4) Remove caster axel and wheel.



5) Using flathead screw drive pop off the end cap on the end of the base of the LiteGait.



7) When loosening the hardware make sure to hold the bolt and square nut to prevent the hardware from falling into the leg of the base.




9) Align the square washer and lock nut with the hole of the tube. Insert the stem bolt into the hole on the base tube.



6) Using 3/4" open ended wrench and 5/16" hex key remove hardware securing caster to base. You may need to utilize a pipe extension to allow for additional leverage when loosening the hardware.




8) Lock each of the new casters prior to installation. Insert stem bolt with lock washer into the new caster fork. Then add the flat washer to the stem bolt of the caster.

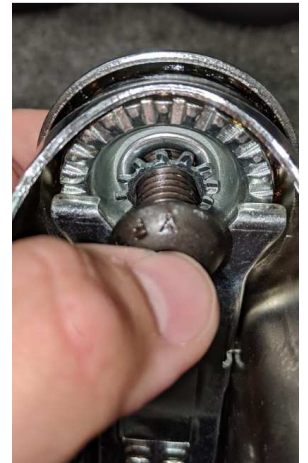

10) Using 3/4" Open ended wrench and 5/16 allen key tighten the hardware until snug to secure the caster to base leg.




11. On the two directional lock casters, tighten the hardware 95% with the caster slightly aligned at a slight angle to the base leg.

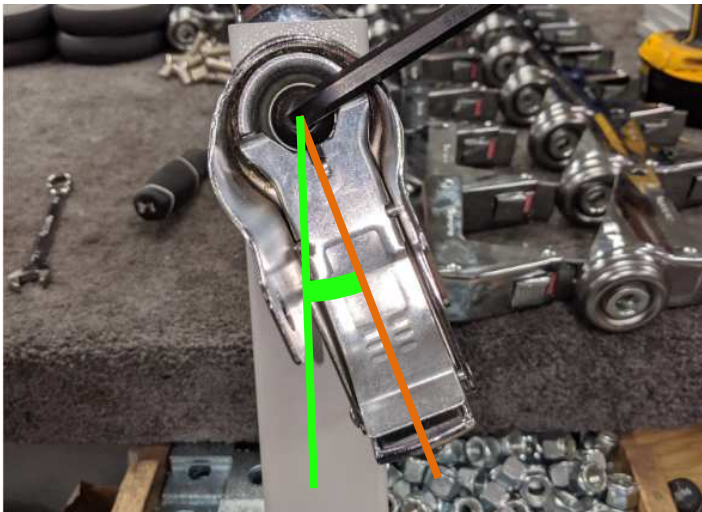

11) Align caster wheel with the hole in the caster fork.

12. Tighten the last 5% to secure the caster inline with the frame and secure the caster to the base leg.

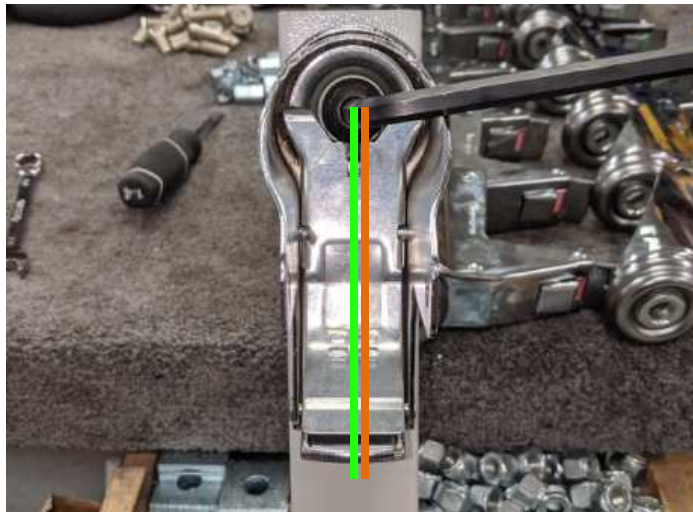

12) Insert carriage bolt into hole with notch cutout. Metal tab on carriage bolt should insert into cutout.



13) Tighten nut until snug using 1/2" open ended wrench.



14) Tighten nut until snug using 1/2" open ended wrench.



15. Reassemble the LiteGait. Please feel free to contact the Service & Parts Department at 1-800-332-9255 Option 1 or email us at Service@LiteGait.com

