




Gaiter Stool



For use with
LiteGait®

Gaiter Stool™ Diagram

1. Adjustable Padded Seat
2. Height Adjustment Lever
3. Adjustable Mounting Arm
4. Stool Base w/Casters
5. Pin Inserts for LiteGait® Base
6. Side Lock

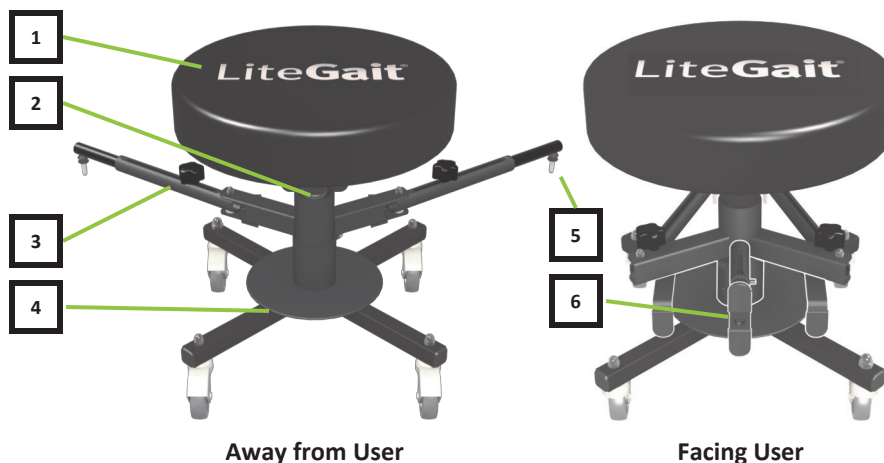


NOTE: Your Gaiter Stool™ may differ from image above

About Your Gaiter Stool™


The Gaiter Stool™ attaches directly to the base of LiteGait® and provides a safe, effective means of assisting patients working over the treadmill or over ground. The Gaiter Stool™ is highly adjustable for improved ease and comfort. The unique design of LiteGait® enables safely supported training of balance, posture and locomotion, while allowing the therapist to observe, evaluate, and correct current gait patterns, manually manipulate limbs, and weight shift to accomplish proper gait patterns. Patients in LiteGait® can easily be moved from over treadmill gait therapy to over ground walking, aiding in the transfer of the newly developed gait patterns.

Adjusting and Attaching Gaiter Stool™


The Gaiter Stool™ can be attached to the LiteGait® via two adjustable arms. The angle of each arm is adjusted by pivoting each arm. When not attached to the LiteGait®, each arm can be collapsed to the stowed position.

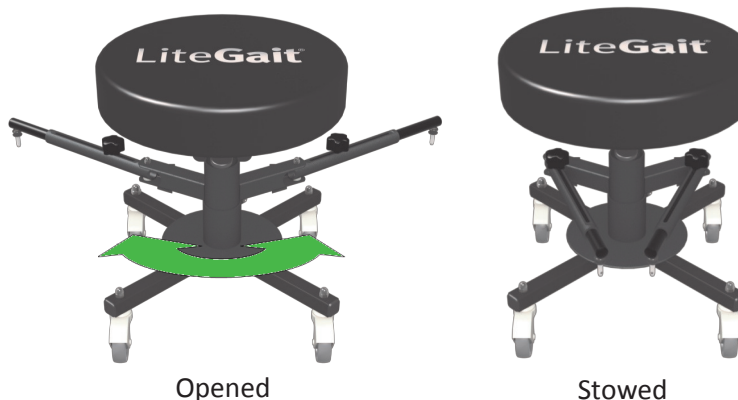

The length of each arm is adjusted by loosening the knob on the top of the arm and pulling out the telescopic arm. Tighten knob to lock arm at designed length.

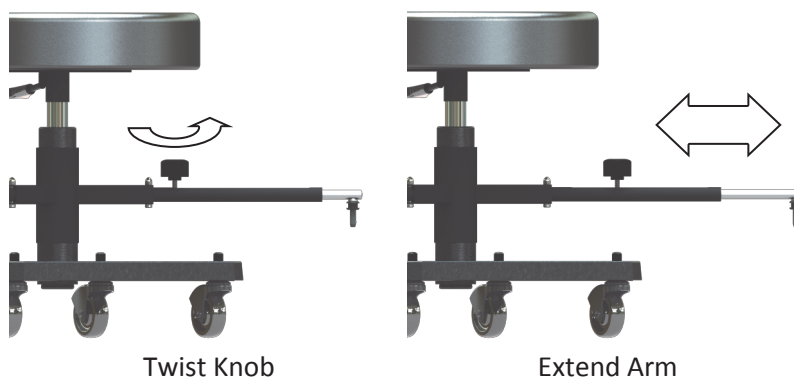

The Gaiter Stool™ can be attached to the LiteGait® via the Gaiter Stool attachment hole on each leg of the LiteGait® base

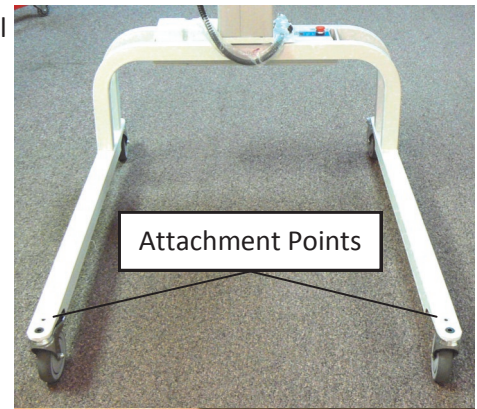


Insert Gaiter Stool™ Pin Insert to attach Gaiter Stool™ to LiteGait® base. To detach Gaiter Stool™, pull up on Pin Insert




The Gaiter Stool™ can be attached to the left or right legs of the base and allowing the therapist to be located on the side of the LiteGait® by using only one of the attachment holes or the side lock option.



Using both attachment holes, the Gaiter Stool™ can be attached to the left and right legs of the base and allowing the therapist to be located behind the user. The angled support of the adjustable arms should open towards the LiteGait®. If the angle support opens away from LiteGait®, then stool may be too close for use.

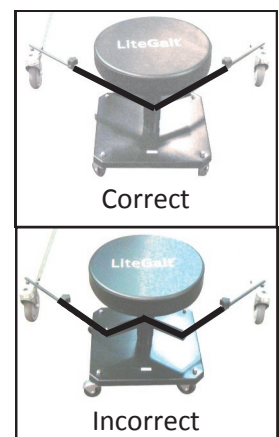

Attaching Gaiter Stool™ using Side Lock feature.

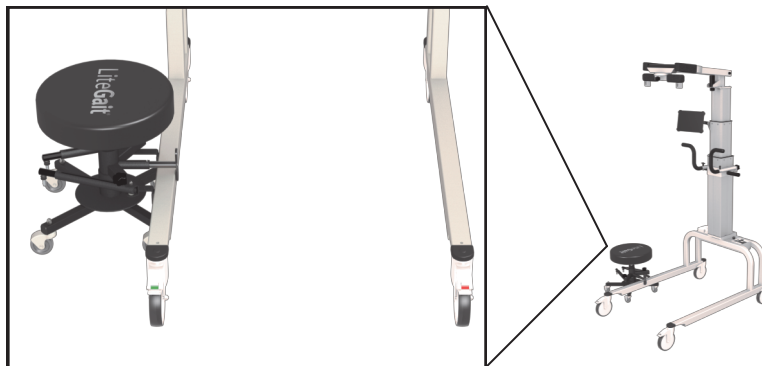

1. Position the Side Lock mechanism towards the outside face of the preferred leg of the LiteGait®.

2. Pull side lock arm over the base leg, away from center stool.

3. Rotate the side lock arm clockwise until vertical.



Unlocked



Locked

Limited Warranty Certificate

Accessories

The Mobility Research warranty covers the following parts only for 1 year against all defects.[†]

- 90 days full replacement warranty
- 1 year on welds and electromechanical parts (excludes wear and tear items, i.e. batteries, soft goods and wheels)

Losses due to work stoppage, lost revenues, damages due to neglect or abuse **ARE NOT** covered by this warranty.
Shipping and handling charges **ARE NOT** covered by this warranty.[†]

THIS WARRANTY AND REMEDIES SET FORTH ABOVE ARE EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES, REMEDIES AND CONDITIONS, WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED. MOBILITY RESEARCH SPECIFICALLY DISCLAIMS ANY AND ALL IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. IF MOBILITY RESEARCH INC. CANNOT LAWFULLY DISCLAIM IMPLIED WARRANTIES UNDER THIS LIMITED WARRANTY, ALL SUCH WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE LIMITED IN DURATION TO THE DURATION OF THIS WARRANTY. NO MOBILITY RESEARCH RESELLER, AGENT, OR EMPLOYEE IS AUTHORIZED TO MAKE ANY MODIFICATION, EXTENSION OR ADDITION TO THIS WARRANTY. MOBILITY RESEARCH IS NOT RESPONSIBLE FOR DIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY BREACH OF WARRANTY OR CONDITION OR UNDER ANY OTHER LEGAL THEORY, INCLUDING BUT NOT LIMITED TO LOST PROFITS, DOWNTIME, GOODWILL, OR DAMAGE TO EQUIPMENT AND PROPERTY.

[†] Warranty excludes damage due to normal wear and tear, tampering with any components, from misuse and abuse, caused by cleaning and acts of God. Warranty does not cover losses due to work stoppage, lost revenue(s), and damages due to neglect. Warranty excludes GaitKeeper Treadmills. Shipping and handling charges are not covered by this warranty. Tampering by persons other than Mobility Research factory authorized personnel is discouraged and will void your warranty.

*In order for us to provide the very best in customer support, please **activate** your warranty by providing the following information. This information will allow us to notify you for product updates, recall information, clinical support, technical support, maintenance information and to receive our E-Newsletter. You may visit our website at <http://www.litegait.com/warranty.html> and submit this form or fill in the information below and mail or fax back. (Keep a copy for your records)*

Model	Serial number	Date of purchase
Facility Name:		Dept. used in
Address		Phone ()
		Fax ()
Clinical Contact Name		email
Maintenance Contact Name		email



Products, Education, and Rehabilitation Solutions

Mobility Research Inc., P. O. Box 3141, Tempe AZ 85280
Phone: 1-800-332-9255 / 480-829-1727 Fax: 480-829-0737
Email: warranty@LiteGait.com

